

REMINGTON®



Paraffin Wax Heat Treatment System for Hands, Feet and Elbows



FEATURES

HS-190 includes:

- Unique Citrus Aromatherapy Paraffin Wax formula
- Full-Size unit designed to accommodate women's shoe size up to 9-1/2
- Soothes and renews hands, elbows and feet
- Moisturizes dry, cracked skin and eases stiff joints

HS-200 includes HS-190 features, plus:

- Insulated hand mitts to extend heat therapy
- Disposable glove liners to protect hand mitt interiors

HS-225 includes HS-190 features, plus:

- Electric massaging hot hand mitts for extended heat and massage therapy
- Disposable glove liners to protect hand mitt interiors

Please Read Prior to Use. Models HS-190, HS-200, HS-225

CONGRATULATIONS...

We at Remington wish to thank you for supporting Remington Products Company, L.L.C. with your purchase. As with all Remington products, your new Paraffin Wax Heat Treatment System has been manufactured with the utmost care and attention to detail and performance to ensure your complete satisfaction.

Relax and refresh as you escape from everyday stresses and strains in a warm bath of aromatic paraffin wax. Just dip your hand, foot or elbow in and out of the wax Heat Treatment System several times to form a "soft" glove of paraffin wax, which releases therapeutic and moisturizing properties deep into the skin and joints. Immediately you can start to feel the heat actually transferring into your skin. Fifteen brief minutes of indulgence are all it takes to refresh, rejuvenate and invigorate your mind and body.

The Remington Paraffin Wax Heat Treatment System provides benefits specifically for your hands, feet and elbows and can help them look and feel better too.

- Soothe and renew hands, feet and elbows
- Moisturize dry and cracked hands, feet and elbows
- Ease stiff joints

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, especially when children are present, basic safety precautions should always be followed, including the following.

READ ALL INSTRUCTIONS BEFORE USING.

DANGER

To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using. See cleaning instructions.
2. Do not use while bathing or in the shower.
3. Do not place or store this appliance where it can fall or be pulled into a tub, sink or dropped into water or other liquid.
4. Do not place or drop into water or other liquid.
5. Do not reach for an appliance that has fallen into water, unplug it immediately.

WARNING

To reduce the risk of burns, electric shock, fire, or injury to persons:

1. This product is for household use only.
2. This product should never be left unattended when plugged in. Always unplug from outlet before adding wax and after using.
3. Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
4. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by Remington.
5. Never operate this product if it has a damaged cord or plug; if it is not working properly; or if it has been dropped or damaged. Return the product to Remington Customer Service for repair or replacement.
6. Keep the cord away from heated surfaces.
7. Never use while sleeping.
8. Never drop or insert any foreign object into the Heat Treatment System chamber.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Do not use an extension cord with this appliance.
11. Do not touch the aluminum side or bottom inside the unit. They are hot.
12. This appliance has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature.
13. It is recommended you use only Remington replacement wax (PW-100). See order form for replacement wax.
14. Do not apply hot wax near eyes or other heat sensitive areas, to any skin injury, eruption, or surface growth.
15. Do not touch hot surfaces.
16. Use this appliance on a level surface.
17. Do not let the cord hang over the edge of the table or counter where it can be pulled on by children or tripped over unintentionally.
18. To disconnect, turn all controls to the off position, then remove plug from the outlet. Allow to cool before handling, cleaning, or storing the appliance.
19. This appliance has no serviceable parts.
20. Do not plug or unplug this product when hands are wet.

21. Do not carry this product by the supply cord or use cord as a handle.
22. Always store this product in a moisture free area.
23. Discontinue use of this appliance if discomfort, irritation, or pain results. Consult your doctor.
24. Do not stand on or in the Heat Treatment System chamber. Use only when seated.
25. Do not pour melted paraffin from the Heat Treatment System chamber.
26. Never use sharp instruments to remove the solidified paraffin from the unit.
27. Do not put water in the Heat Treatment System chamber or immerse the unit in water.
28. Do not allow the level of paraffin in the Heat Treatment System chamber to fall below the "min" line, or above the "max" line. With the unit filled to the "max" line, it takes approximately 8–10 uses (5 dips per hand per use) for the wax to fall below the "min" line.
29. Do not heat paraffin on a range, stove, or in a microwave, or in plastic packaging material.
30. Do not use any additives in the paraffin.


PLEASE NOTE: This appliance should not be used in the presence of open cuts or wounds, inflammatory skin conditions, growth peripheral vascular disease where circulation is impaired, acute inflammation, or where sensation of the extremity is reduced or absent. If there is any question about peripheral vascular disease or decreased sensation of the extremities, consult a physician or physical therapist before using the Paraffin Wax Heat Treatment System. Discontinue use if dermatitis due to paraffin sensitivity occurs. Discontinue use if wax feels too hot or cool, which could indicate health problems with the user.

SAVE THESE INSTRUCTIONS

OPERATING INSTRUCTIONS

How to Fill the Spa

1. Remove the lid. Set all documents aside.
2. Place the unit on a hard, flat, heat-resistant surface.
3. Remove all plastic packaging material from the paraffin wax. Break all 4 lbs. of wax into small pieces, then place in Heat Treatment System chamber. The smaller the pieces, the faster the heat up time.
4. Make sure unit is in the "off" position. Plug unit into an electrical outlet of proper voltage.
5. This product has a polarized plug (one blade is wider than the other). As a safety feature, this plug will fit in a polarized outlet only one way. If the plug does not fully fit in the outlet, reverse the plug. If it still does not fit in the outlet, contact a qualified electrician to install a polarized outlet. Never try to defeat this safety feature.
6. Replace the lid to retain heat. Turn on the unit by sliding the switch button to the "on" position and wait until the wax is visibly melted (approximately 90 minutes). Once wax is melted, do not move unit as spillage may occur!
7. With the unit filled to the "max" line, it takes approximately 8 – 10 uses (5 dips per hand per use) for the wax to fall below the "min" line.

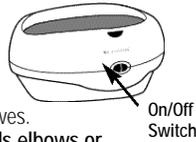

How to Dip Hands

1. Wash and thoroughly dry hands, remove jewelry and roll up sleeves.
For ultimate softness, add your favorite lotion to your hands, elbows, or feet and rub in gently prior to dipping. Important: Do not use a heat activated or heat sensitive lotion. NOTE: Lid should be on the unit during heat up time. Lid should be placed back on the unit after each use.
2. Test the temperature of the wax by dabbing a small amount on the inside of your wrist. If wax seems too hot, simply unplug the unit and allow wax to cool slightly before dipping.
3. With fingers spread slightly apart, slowly immerse your hand into the Heat Treatment System chamber. Then immediately, but gently, remove your hand from the Heat Treatment System chamber. Let the paraffin solidify on your hand. Wait at least three seconds (until the paraffin no longer appears shiny) before redipping. **Be careful to avoid touching the sides or the bottom of the Heat Treatment System chamber.**
4. Continue to dip your hand in and out of the wax, following Step 3, to build up layers of paraffin. For the best results, do this at least five times. The more layers, the longer the heat will last.
5. If desired, cover your hands with the plastic glove liners (HS-200 and HS-225 only) and place your hands in the insulated (HS-200 only) or massaging (HS-225 only) heated mitts. This retains heat for better results. Do not use massaging hot mitts (HS-225 only) for more than 15 minutes.
6. After massage is completed, remove the hand mitts from your hands. Remove the plastic glove liners and discard. Peel the hardened wax from your hands. **It is recommended that the hardened wax be discarded due to the buildup of dead skin cells**

attached to the wax. However, you may put the hardened wax back in the chamber, if desired.

7. Massage in the oils left on your hand after the wax is peeled off to get the full invigorating and rejuvenating benefits of your Paraffin Wax Heat Treatment System.
8. As soon as you have completed your wax treatment, turn the unit off and unplug it. If you discard the used wax after each use, we recommend you clean out the unit completely after 40 uses. If you re-use the wax after every use, we recommend you clean out the unit completely after 20 uses. (See "Cleaning Instructions").
9. If, at this time, cleaning is not necessary simply allow wax to harden completely and then cover unit and put away until next usage. When you are ready to use the Wax Heat Treatment System again, simply remove from storage, plug unit into a polarized position and wait until wax is visibly melted (approximately 90 minutes). See "Operating Instructions" for further instructions.

To Treat Feet or Elbows

Follow steps 1-4 as outlined for hands. Glove liners and mitts are intended for use with hands only.

NOTE: Allow the paraffin to remain in place for 10-15 minutes before removing, longer if more layers exist or if the insulated mitts (HS-200 only) were used.

OPERATING INSTRUCTIONS FOR MASSAGING HOT MITTS (HM-300)

Refer to instruction book included with Salon Style Hot Mitts.

CLEANING INSTRUCTIONS FOR PARAFFIN WAX HEAT TREATMENT SYSTEM

1. If you discard the used wax after each use, we recommend you clean out the unit completely every 40 uses. If you re-use the wax after every use, we recommend you clean out the unit completely every 20 uses.
2. Turn unit off.
3. **Before unplugging the Wax Heat Treatment System**, allow paraffin to remain undisturbed for a few minutes.
4. Unplug unit and let the wax "harden" overnight.
5. Right before cleaning, plug unit in and turn the unit on again for approximately 1-2 minutes, then turn unit off and unplug from the electrical outlet.
6. Remove hardened wax block by using several paper towels and slowly wiping along the bottom of the Heat Treatment System chamber, collecting any sediment in the paper towels. Slowly bring paper towels out of the paraffin, being careful not to let the sediment drop back in. A non-abrasive all-purpose cleaner may be used to clean the inside of the chamber.
7. Do not use abrasives when cleaning the inside of the Heat Treatment System chamber as it will damage the finish.
8. To clean the outside of the unit, remove any solidified pieces of paraffin and use a non-abrasive all-purpose cleaner.
9. Do not use sharp objects to remove solidified paraffin.
10. Do not scratch or damage Heat Treatment System chamber.

For Replacement Paraffin Wax or Glove Liners, see order form or call 1-800-736-4648 in U.S. or 1-800-268-6531 in Canada.

FULL TWO YEAR WARRANTY

Remington Products Company, L.L.C. warrants this product against any defects that are due to faulty material or workmanship for a two-year period from the original date of consumer purchase. This warranty does not include damage to the product resulting from accident or misuse.

If the product should become defective within the warranty period, we will repair it or elect to replace it free of charge. We will return your product, transportation charges prepaid, provided it is delivered prepaid to: Remington Products Company, L.L.C., 60 Main Street, Bridgeport, CT 06604, Attention: Customer Service Department. In Canada: Remington Products (Canada) Inc., 455 Cochrane Drive, Unit #24, Markham, Ontario L3R 9R4, Attention: Customer Service. Include a letter stating the reason for the return.

REMINGTON SHALL NOT BE RESPONSIBLE FOR ANY INCIDENTAL, SPECIAL OR CONSEQUENTIAL DAMAGES.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state. Some states do not allow the exclusion or limitation of incidental, special or consequential damages.

REMINGTON MONEY BACK OFFER

If within 90 days after you have purchased your Remington Paraffin Wax Heat Treatment System, you are not satisfied and would like a refund, take back the item and sales slip indicating purchase price and date of purchase to the retailer from whom the item was purchased. Remington will reimburse all retailers who accept the item within 90 days from the date of purchase. If you have any questions concerning the money back guarantee, please call 1-800-736-4648 in the U.S. only or 1-905-470-9400 in Canada.



This product will bear one or more compliance marks, Listed by Underwriters Laboratories, Inc. 120 VAC, 60 Hz



Questions or Comments: Please call 1-800-736-4648 in U.S. only or 1-905-470-9400 in Canada.



Visit us at <http://www.remington-products.com>

Remington. Better products. Better value. Guaranteed.

REMINGTON

PRODUCTS COMPANY, L.L.C.

60 Main Street Bridgeport, CT 06604

Remington Products (Canada) Inc., 455 Cochrane Drive,
Unit #24, Markham, Ontario L3R 9R4.

Remington and Spa Therapy Collection design are
trademarks of Remington Corporation, L.L.C.
©1999 R.P.C., L.L.C. Rev. 9/99 Job#99-092 Part# 13084